



# Pacific Medical Catalogue

TABLE OF CONTENTS

**Electrosurgery** ..... 2

    ESU Devices ..... 2

    ESU Device Mounts & Accessories ..... 8

    Handpieces ..... 11

    Sheaths ..... 15

    Return Plates & Scratch Pads ..... 17

    Disposable Electrodes ..... 20

    Reusable Electrodes ..... 31

    Bipolar Forceps & Accessories ..... 37

    Smoke Evacuators & Accessories ..... 43

    Miscellaneous & Spare Parts ..... 51

    Cautery ..... 53

        Power Handles & Replacement Tips ..... 53

        High Temperature Cauteries ..... 58


        Low Temperature Cauteries ..... 62

        Sleeves ..... 64

    Ophthalmic ..... 65

# ELECTROSURGERY

## ESU DEVICES

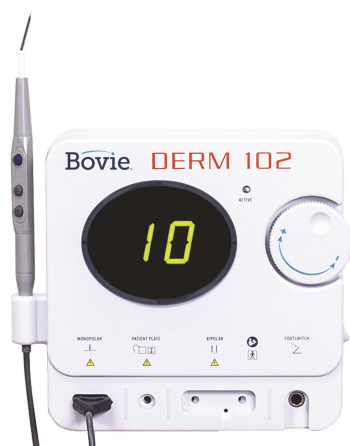

### BOVIE DERM 102 HIGH FREQUENCY DESICCATOR

Most procedures performed by a High Frequency Desiccator utilize ten watts of power or less. The DERM 102 is an economical, 10 watt high frequency desiccator, thoughtfully designed to give practitioners the ability to perform in-office, minor skin procedures such as the removal of skin tags, benign lesions and premalignant lesions. The Bovie Derm 102 unit allows practitioners to provide better patient care in the office, at a sensible price, saving the patient time and money, while increasing their practice revenue.

The Bovie DERM 102 High Frequency Desiccator offers ten watts of monopolar power, adjustable in 1/10th watt increments to desiccate and fulgurate tissue. The Derm 102 unit power can be adjusted using the 3-Button Power Control Handpiece or by using the rotary dial on the unit. The unit features a large, bright, LED display allowing excellent visibility across the procedure room.

For the practitioner who wants to also use bi-polar electrosurgery along with high frequency desiccation, the Bovie DERM 102 High Frequency Desiccator provides 10 watts of bipolar power, adjustable in 1/10th watt increments. To utilize this feature the bipolar cable, a bipolar forcep and footswitch activation (all sold separately) are required.

The Bovie DERM 102 is packaged ready to use with a total of 18 disposable dermal tips

#### Features:

- Large illuminated display
- Audible tones
- Monopolar & bipolar outputs
- Higher resolution in bipolar mode
- Rotary power dial
- Isolated circuitry

#### Specification:

##### Output Power

To 10 watts (0-10 watts in 1/10th watt increments on both monopolar & bipolar settings)

##### Dimensions

Height: 20.8cm  
Width: 18.4cm  
Depth: 6.9cm  
Weight: 1.4kg

**Main Current**

320mA Max

**Line Voltage**

100-240VAC

**Line Frequency**

50-60Hz

**Warranty**

Industry leading 4 year warranty

**Kit Inclusions**

Derm 102 high frequency desiccator x 1  
A902 power control pencil x 1  
A910 disposable sheath x 2  
A804 non-sterile single use sharp dermal tips x 5  
A805 sterile single use sharp dermal tips x 2  
A806 non-sterile single use blunt dermal tips x 5  
A807 sterile single use blunt dermal tips x 2  
A806DE non-sterile single use DERM-Elite premium blunt dermal tips x 2  
A807DE sterile single use DERM-Elite premium blunt dermal tips x 2  
A837 wall mount kit x 1  
AC power cord x 1  
User guide x 1

**ORDER CODES:**

**BOVDERM102**

**MAINTENANCE:** Please make sure your device is checked and serviced yearly by an authorised technician.

[Click For More Information](#)

## BOVIE DERM 942 HIGH FREQUENCY DESICCATOR

The Bovie DERM 942 High Frequency Desiccator and its patent pending handpiece are designed and engineered to be reliable and durable. Power is adjustable in 1/10th watt increments in the first 10 watts, in monopolar and bipolar. The Bovie DERM 942 is a ground reference unit, which means no grounding pad is required.

The Bovie Derm 942 High Frequency Desiccator and its patent pending handpiece are designed and engineered to be reliable and durable. Safety features such as digital error detection continuously monitors every aspect of the units output and at any sign of problem instantly disables output and displays an error code. The Bovie Derm 942 unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal



heat. Post-procedure and long-term clean-up are simplified by its non-absorbing, non-staining, high impact plastic case. Power is adjustable in 1/10th watt increments in the first 10 watts, in monopolar and bipolar. The Bovie DERM 942 is a ground reference unit, which means no grounding pad is required reducing overall device ownership cost.

**Features:**

Large illuminated display  
Self test circuits  
Audible tones  
Monopolar & bipolar outputs  
Higher resolution in bipolar mode  
Rotary power dial  
Isolated circuitry  
Digital error detection

**Specification:****Output Power**

To 40 watts (0-10 watts in 1/10th watt increments on both monopolar & bipolar settings)

**Dimensions**

Height: 18.8cm  
Width: 22.8cm  
Depth: 10.5cm  
Weight: 2.26kg

**Output Frequency**

328-492kHz

**Line Voltage**

100-240VAC

**Line Frequency**

50-60Hz

**Warranty**

Industry leading 4 year warranty

**Kit Inclusions**

A942 electrosurgical generator x 1  
A902 power control pencil x 1  
A910 disposable sheath x 2  
A804 non-sterile single use sharp dermal tips x 5  
A805 sterile single use sharp dermal tips x 2  
A806 non-sterile single use blunt dermal tips x 5  
A807 sterile single use blunt dermal tips x 2  
A806DE non-sterile single use DERM-Elite premium blunt dermal tips x 2  
A807DE sterile single use DERM-Elite premium blunt dermal tips x 2  
A837 wall mount kit x 1  
AC power cord x 1  
User guide x 1

**ORDER CODES:**

**BOVA942**

**MAINTENANCE: Please make sure your device is checked and serviced yearly by an authorised technician.**

[Click For More Information](#)



## BOVIE BANTAM PRO HIGH FREQUENCY DESICCATOR

The Bovie Bantam PRO is a high frequency desiccator that also has cutting ability. The unit provides 50 watts of cutting power and coagulation modes, plenty of power for most procedures performed in the doctor's office. Power is adjustable in 1/10th watt increments from 0 to 10 watts, for the most accurate output available.

The Bovie Bantam PRO (A952) is a High Frequency Desiccator with cutting ability. You might call it a split personality. One half of the unit is a high frequency desiccator, an upgrade to our very popular Aaron 940. The other half is a Bovie Generator, giving you 50 watts of cutting power and coagulation capabilities. The Bovie Bantam PRO gives you plenty of power for most cutting procedures performed in the doctor's office. The unit comes with a suite of accessories and is ready to use straight from the box.

### Features:

- Large illuminated display
- Self test circuits
- Audible tones
- Monopolar & bipolar outputs
- Higher resolution in bipolar mode
- Rotary power dial
- Isolated circuitry
- Digital error detection

### Specification:

#### Output Power

Cut: 50 watts (0-10 watts in 1/10th watt increments)  
 Blend: 50 watts  
 Coagulation: 50 watts  
 Fulguration: 40 watts  
 Bipolar & Micro Bipolar: 40 watts

#### Dimensions

Height: 18.8cm  
 Width: 22.8cm  
 Depth: 10.5cm  
 Weight: 2.26kg

#### Output Frequency

304-451kHz

#### Line Voltage

100-240VAC

#### Line Frequency

50-60Hz

### Warranty

Industry leading 4 year warranty

**Kit Inclusions**

A952 electrosurgical generator x 1  
A902 power control pencil x 1  
A910 disposable sheath x 2  
A910ST sterile disposable sheath x 2  
A804 non-sterile single use sharp dermal tips x 5  
A805 sterile single use sharp dermal tips x 2  
A806 non-sterile single use blunt dermal tips x 5  
A807 sterile single use blunt dermal tips x 2  
A806DE non-sterile single use DERM-Elite premium blunt dermal tips x 2  
A807DE sterile single use DERM-Elite premium blunt dermal tips x 2  
ES01 sterile single use blade electrode x 3  
ES02 sterile single use needle electrode x 1  
ES20 sterile single use ball electrode x 1  
ESRE single use split grounding pad x 5  
A1252C reusable grounding cord x 1  
A837 wall mount kit x 1  
AC power cord x 1  
User guide x 1

**ORDER CODES:  
BOVA952**

**MAINTENANCE: Please make sure your device is checked and serviced yearly by an authorised technician.**

[Click For More Information](#)

## BOVIE SPECIALIST PRO HIGH FREQUENCY ELECTROSURGICAL GENERATOR

The Bovie Specialist PRO is the affordable solution for the busy physician office and specialty facility. The multi-purpose unit features a modern push button design and is capable of most monopolar and bipolar surgical procedures. The Specialist PRO offers 5 energy outputs with up to 120 watts of power; modalities include Cut, Blend, Coagulation, Fulguration and Standard Bipolar. The unit is designed with dependable safety features like BovieNEM™ (Neutral Electrode Monitoring) pad sensing technology, BovieDED™ (Digital Error Detection), audible alert tones and self-test circuits. The Bovie Specialist | PRO satisfies the physician's electrosurgical needs with safety, reliability and convenience.

The Bovie Specialist PRO is the affordable solution for the busy physician office or specialty facility. Capable of most monopolar and bipolar surgical procedures and packed with dependable safety features such as digital error detection & isolated circuitry the Bovie Specialist PRO satisfies the physician's surgical demands with safety, reliability and convenience. When paired with the Bovie Smoke Shark II Smoke Evacuator this unit is a total system solution for procedures and facilities that require high performance with smoke evacuation.



**Features:**

Large illuminated display  
Ample raised push button power-controls providing rapid power adjustments  
Self test circuits  
Audible tones  
Monopolar & bipolar outputs  
Discreet Outputs FCFS (First Come First Served)  
Isolated circuitry  
Digital error detection

**Specification:****Output Power**

Cut: 120 Watts  
Blend: 90 Watts  
Coagulation: 80 Watts  
Fulguration: 40 Watts  
Bipolar: 30 Watts

**Dimensions**

Height: 15.2cm  
Width: 26cm  
Depth: 30.5cm  
Weight: 6.5kg

**Output Frequency**

350-600kHz

**Line Voltage**

100-240VAC

**Line Frequency**

50-60Hz

**Warranty**

Industry leading 4 year warranty

**Kit Inclusions**

A1250S electrosurgical generator x 1  
ESP1 sterile disposable pencil x 1  
ES01 sterile single use blade electrode x 1  
ES02 sterile single use needle electrode x 1  
ES20 sterile single use ball electrode x 1  
ESRE single use split grounding pads x 5  
A1252C reusable grounding cord x 1  
AC power cord x 1  
User guide x 1

**ORDER CODES:**


**BOVA1250S**

**MAINTENANCE: Please make sure your device is checked and serviced yearly by an authorised technician.**

[Click For More Information](#)

# ESU DEVICE MOUNTS & ACCESSORIES

## BOVIE WALL MOUNT KIT

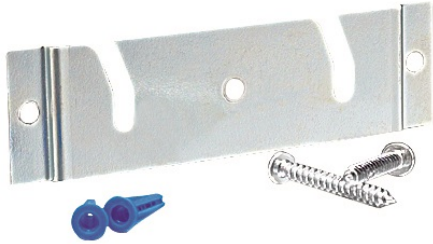

A mount to allow easy wall storage of Bovie devices.

### Features:

Includes mounting plate & hardware  
Suitable for Bantam PRO, Derm 942 & Derm 101 units

**ORDER CODES:**  
**BOVA837**

[Click For More Information](#)

## BOVIE MOBILE STAND

A mobile stand to allow easy transfer of Bovie devices.

### Features:

Includes mounting plate.  
Suitable for Bantam PRO, Derm 942 & Derm 101 units

**ORDER CODES:**  
**BOVA812**

[Click For More Information](#)





## **BOVIE MOBILE CART FOR ESU & SMOKE EVACUATOR**

A three tier mobile cart to suit Bovie ESU & Smoke Shark smoke evacuator,

### **Features:**

Includes three tiers  
Smooth rolling lockable castors  
Suitable for Bovie Specialist PRO & Smoke Shark 2.

### **ORDER CODES: BOVESMS2**

[Click For More Information](#)

## **BOVIE FOOT SWITCH**

A Bovie foot switch to suit the Derm 102, Derm 942 & Bantam Pro units.

### **ORDER CODES: BOVA803**

[Click For More Information](#)



## BOVIE MONOPOLAR / BIPOLAR FOOT SWITCH



Operates both monopolar & bipolar functions on the A1250S Specialist Pro ESU.

**ORDER CODES:**  
**BOVBV1253B**

[Click For More Information](#)

## BOVIE POWER CABLE

A power cable to suit all Bovie units.

**ORDER CODES:**  
**BOVPC**

\*\* Please note power cord and plug colour may vary \*\*

[Click For More Information](#)



# HANDPIECES

## BOVIE DISPOSABLE PUSH BUTTON PENCIL WITH HOLSTER & BLADE BX40



An disposable push button pencil including a holster and standard blade electrode.

**ORDER CODES:**  
**BOVESP1H**

[Click For More Information](#)

## AUTOCLAVABLE THREE BUTTON PENCIL NEW

An Autoclavable three button pencil to suit Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952) Bovie generators.

**ORDER CODES:**  
**BOVA902**



[Click For More Information](#)

## AUTOCLAVABLE THREE BUTTON PENCIL OLD



An Autoclavable three button pencil to suit A900 & A940 Bovie generators

**ORDER CODES:**  
**BOVA901**

[Click For More Information](#)

## BOVIE AUTOCLAVABLE REUSABLE PENCIL

An Autoclavable reusable pencil validated to 100 autoclave cycles.

**ORDER CODES:**  
**BOVESPR2**



[Click For More Information](#)

## **BOVIE SUCTION COAGULATOR DISPOSABLE HANDSWITCH 12FR 3M CABLE PK10**



Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

This coagulator is 12 french gauge in size, operated by handswitch and has a 3M cable.

**ORDER CODES:  
BOV SCH12**

[Click For More Information](#)

## **BOVIE SUCTION COAGULATOR DISPOSABLE FOOTSWITCH 10FR 3M CABLE PK10**

Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

This coagulator is 10 french gauge in size, operated by footswitch and has a 3M cable.

**ORDER CODES:  
BOV SCF10**

[Click For More Information](#)



## **BOVIE SUCTION COAGULATOR DISPOSABLE HANDSWITCH 8FR 3M CABLE PK10**



Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

This coagulator is 8 french gauge in size, operated by handswitch and has a 3M cable.

**ORDER CODES:  
BOVSCH08**

[Click For More Information](#)

# SHEATHS

## BOVIE DISPOSABLE HANDPIECE SHEATH NON-STERILE BX100



Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie generators, as well as most competitor generators.

Bovie handpiece sheath are used in conjunction with electrosurgery pencils to provide a cost-effective alternative to continuous sterilization. The use of handpiece sheaths significantly increases the life of the reusable pencils.

**ORDER CODES:**  
**BOVA910**

[Click For More Information](#)

## BOVIE DISPOSABLE HANDPIECE SHEATH STERILE BX25

Bovie handpiece sheath are used in conjunction with electrosurgery pencils to provide a cost-effective alternative to continuous sterilization. The use of handpiece sheaths significantly increases the life of the reusable pencils.

Bovie handpiece sheath are used in conjunction with electrosurgery pencils to provide a cost-effective alternative to continuous sterilization. The use of handpiece sheaths significantly increases the life of the reusable pencils.

**ORDER CODES:**  
**BOVA910ST**

[Click For More Information](#)





# RETURN PLATES & SCRATCH PADS



## BOVIE REUSABLE METAL GROUNDING PLATE WITH CABLE

All Bovie patient return electrodes are manufactured to the highest standards and come individually packaged.

Reusable metal grounding plate is used to safely remove electrosurgical current from the patient. Bovie recommends always using a grounding plate for maximum patient safety.

This product is suitable for use with the Derm 101, Derm 102 & Derm 942 (A942) units.

**ORDER CODES:**  
**BOVA802EU**

[Click For More Information](#)

## BOVIE DISPOSABLE SPLIT ADULT RETURN ELECTRODE WITH 3M CABLE BX50

All Bovie patient return electrodes are manufactured to the highest standards and come individually packaged.

Bovie disposable split adult return electrodes are used to safely remove electrosurgical current from the patient and are manufactured with an adhesive hydrogel for secure patient application. Bovie recommends always using a split pad for maximum patient safety.

**ORDER CODES:**  
**BOVESREC**

[Click For More Information](#)



## BOVIE DISPOSABLE SPLIT ADULT RETURN ELECTRODE BX50



All Bovie patient return electrodes are manufactured to the highest standards and come individually packaged.

Bovie disposable split adult return electrodes are used to safely remove electrosurgical current from the patient and are manufactured with an adhesive hydrogel for secure patient application. Bovie recommends always using a split pad for maximum patient safety.

**ORDER CODES:**  
**BOVESRE**

[Click For More Information](#)

## BOVIE REPLACEMENT METAL GROUNDING PLATE

All Bovie patient return electrodes are manufactured to the highest standards and come individually packaged.

Reusable metal grounding plate is used to safely remove electrosurgical current from the patient. Bovie recommends always using a grounding plate for maximum patient safety.

This product is suitable for use with the A1204C & 1254C cables.

**ORDER CODES:**  
**BOVA1204P**



[Click For More Information](#)



## BOVIE REUSABLE GROUNDING CABLE

All Bovie patient electrode cables are manufactured to the highest standards and come individually packaged.

Reusable grounding cable is used to connect to a return patient electrodes. Bovie recommends always using a grounding plate for maximum patient safety.

This product is suitable for use with ESRE & ESRS return electrodes. It is also suitable for use with Specialist Pro (A1250S) & Derm 942 (A942) electrosurgical units.

**ORDER CODES:**  
**BOVA1252C**


[Click For More Information](#)

## BOVIE SCRATCH PADS BX40

All Bovie patient electrode cables are manufactured to the highest standards and come individually packaged.

Bovie Scratch pads are used to clean stubborn build up from Electrodes. They come in a box of 40 units.

**ORDER CODES:**  
**BOVESSP**

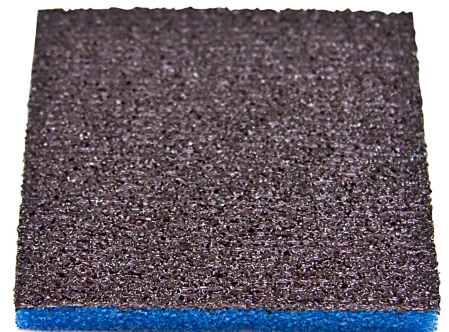

[Click For More Information](#)

# DISPOSABLE ELECTRODES



## BOVIE DISPOSABLE DERMAL ELECTRODE TIPS SHARP NON-STERILE BX100

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVA804**

[Click For More Information](#)

## BOVIE DISPOSABLE DERMAL ELECTRODE TIPS SHARP STERILE BX50

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVA805**

[Click For More Information](#)





## **BOVIE DISPOSABLE DERMAL ELECTRODE TIPS BLUNT ANGLED NON-STERILE BX100**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVA806**

[Click For More Information](#)

## **BOVIE DISPOSABLE DERMAL ELECTRODE TIPS BLUNT ANGLED STERILE BX50**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVA807**

[Click For More Information](#)





### **BOVIE DISPOSABLE DERM-ELITE ROUND ELECTRODE TIPS BLUNT ANGLED NON-STERILE BX100**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These round Derm-Elite electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVA806DE**

[Click For More Information](#)

### **BOVIE DISPOSABLE DERM-ELITE ROUND ELECTRODE TIPS BLUNT ANGLED STERILE BX50**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These round Derm-Elite electrodes are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVA807DE**

[Click For More Information](#)



## **BOVIE DISPOSABLE STANDARD BLADE ELECTRODE BX50**



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable blade electrode are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVES01**

[Click For More Information](#)

## **BOVIE DISPOSABLE RESTICK II COATED BLADE ELECTRODES BX12**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable Restick II coated blade electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVES01T**

[Click For More Information](#)



## BOVIE DISPOSABLE EXTENDED BLADE ELECTRODES BX25



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 15.24cm disposable extended blade electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES04**

[Click For More Information](#)

## BOVIE DISPOSABLE STANDARD NEEDLE ELECTRODES BX25

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES02**

[Click For More Information](#)



## BOVIE DISPOSABLE RESTICK II COATED NEEDLE ELECTRODES BX12



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable Restick II coated needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES02T**

[Click For More Information](#)

## BOVIE DISPOSABLE MICRO DISSECTION NEEDLE ELECTRODE BX5

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 4.5cm disposable micro dissection needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES61**

[Click For More Information](#)



## **BOVIE DISPOSABLE RESTICK II COATED MODIFIED NEEDLE ELECTRODES BX12**



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable Restick II coated modified needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVES38T**

[Click For More Information](#)

## **BOVIE DISPOSABLE SUPER CUT MODIFIED FINE NEEDLE ELECTRODES BX5**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 4.5cm disposable super cur modified fine needle electrode are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVES63**

[Click For More Information](#)



## **BOVIE DISPOSABLE IMA / ENT STANDARD BLADE ELECTRODES BX25**



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable IMA / ENT standard blade electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES37**

[Click For More Information](#)

## **BOVIE DISPOSABLE IMA / ENT STANDARD NEEDLE ELECTRODES BX25**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable IMA / ENT standard needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES38**

[Click For More Information](#)



## **BOVIE DISPOSABLE IMA / ENT EXTENDED NEEDLE ELECTRODES BX25**



All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These 6.99cm disposable IMA / ENT extended needle electrodes are suitable for use with all Bovie units including the Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVES40**

[Click For More Information](#)

## **BOVIE DISPOSABLE H-TYPE ELECTRODE TIPS SHARP ANGLED NON-STERILE BX100**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These H-Type electrodes sharp non-sterile are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:**  
**BOVH10012**

[Click For More Information](#)





## **BOVIE DISPOSABLE H-TYPE ELECTRODE TIPS SHARP STERILE BX50**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These H-Type electrodes sharp sterile are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVH10008**

[Click For More Information](#)

## **BOVIE DISPOSABLE H-TYPE ELECTRODE TIPS BLUNT ANGLED NON-STERILE BX100**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These H-Type electrodes blunt non-sterile are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVH10112**

[Click For More Information](#)





## **BOVIE DISPOSABLE H-TYPE ELECTRODE TIPS BLUNT ANGLED STERILE BX50**

All electrodes utilize the standard 2.3mm stainless steel shaft electrodes are manufactured to the highest standards and come individually packaged.

These H-Type electrodes blunt sterile are suitable for use with the Bovie Derm 101, Derm 102, Derm 942 (A942) & Bantam PRO (A952). They are also suitable for use with any pencil that accepts the standard 2.3mm stainless steel shaft.

**ORDER CODES:  
BOVH10108**

[Click For More Information](#)

# REUSABLE ELECTRODES

## BOVIE REUSABLE EPILATION NEEDLE ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA834**

[Click For More Information](#)

## BOVIE REUSABLE STRAIGHT NEEDLE ELECTRODE

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable straight needle electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA833**

[Click For More Information](#)



## BOVIE REUSABLE STANDARD NEEDLE ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable angled needle electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVES02R**

[Click For More Information](#)

## BOVIE REUSABLE 5MM EXTENDED BALL ELECTRODE

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.



Reusable 5mm extended ball electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVES07R**

[Click For More Information](#)

## BOVIE REUSABLE ANGLED BALL ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable angled ball electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA831**

[Click For More Information](#)

## BOVIE REUSABLE STANDARD BLADE ELECTRODE

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable standard blade electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVES01R**

[Click For More Information](#)



## BOVIE REUSABLE STANDARD SEGAL ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable standard segal electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA811**

[Click For More Information](#)

## BOVIE REUSABLE 1/4" TUNGSTEN LOOP ELECTRODE

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.



Reusable 1/4" loop electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVES22R**

[Click For More Information](#)

## BOVIE REUSABLE ANGLED NEEDLE ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable angled needle electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA830**

[Click For More Information](#)

## BOVIE REUSABLE 3/8" TUNGSTEN LOOP ELECTRODE

All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.



Reusable 3/8" loop electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVES24R**

[Click For More Information](#)

## BOVIE REUSABLE 3/32" BALL ELECTRODE



All reusable electrodes utilize the standard 2.3mm stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil.

Reusable 3/32" ball electrode utilize the standard 2.3mm stainless steel shaft. The moulded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for 5 - 25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

**ORDER CODES:**  
**BOVA832M**

[Click For More Information](#)

# BIPOLAR FORCEPS & ACCESSORIES



## BOVIE REUSABLE BIPOLAR FORCEPS CORD FOR NEW UNITS

This cable allows connection of current generation Bovie electrosurgical units to bipolar forceps.

This cable allows connection of current Bovie electrosurgical units to bipolar forceps. The cable is 3M long and suited to devices such as the Derm 101, Derm 102, Derm 942 (A942), Bantam Pro (A952) & Specialist Pro (A1250S)

**ORDER CODES:**  
**BOVA827V**

[Click For More Information](#)

## BOVIE 184MM FINE TIP BIPOLAR FORCEP

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Bovie 184mm fine tip bipolar forcep are guaranteed for twenty autoclave cycles. Bipolar forceps require an electrosurgery footswitch for activation, see all Bovie accessories for the available range. For care and handling, refer to the Instructions for Use in the package insert.

### Features:

- High quality materials
- Guaranteed for 20 autoclave cycles
- Require an electrosurgery footswitch for activation

**ORDER CODES:**  
**BOVA820**

[Click For More Information](#)



## BOVIE 177MM STRAIGHT FINE TIP BIPOLAR FORCEP



Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Bovie 177mm straight fine tip bipolar forcep are guaranteed for twenty autoclave cycles. Bipolar forceps require an electrosurgery footswitch for activation, see all Bovie accessories for the available range. For care and handling, refer to the Instructions for Use in the package insert.

### Features:

- High quality materials
- Guaranteed for 20 autoclave cycles
- Require an electrosurgery footswitch for activation

**ORDER CODES:**  
**BOVA822**

[Click For More Information](#)

## BOVIE 114MM STRAIGHT MICRO TIP BIPOLAR FORCEP

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Bovie 114mm straight micro tip bipolar forcep are guaranteed for twenty autoclave cycles. Bipolar forceps require an electrosurgery footswitch for activation, see all Bovie accessories for the available range. For care and handling, refer to the Instructions for Use in the package insert.

### Features:

- High quality materials
- Guaranteed for 20 autoclave cycles
- Require an electrosurgery footswitch for activation

**ORDER CODES:**  
**BOVA826**

[Click For More Information](#)



## **OLSEN 178MM CUSHING SMOOTH STRAIGHT 1.5MM TIPS INSULATED BIPOLAR FORCEP**



Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Olsen 178mm cushing smooth straight 1.5mm tips are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### **Features:**

High quality materials  
Insulated  
Suitable for use with all electrosurgical devices

**ORDER CODES:**  
**BOV10-11621**

[Click For More Information](#)

## **OLSEN 127MM ADSON SMOOTH 1.0MM TIPS INSULATED BIPOLAR FORCEP**

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Olsen 127mm adson smooth 1.0mm tips are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### **Features:**

High quality materials  
Insulated  
Suitable for use with all electrosurgical devices

**ORDER CODES:**  
**BOV10-10821**

[Click For More Information](#)





## BOVIE REUSABLE BIPOLAR FORCEPS CORD FOR OLD UNITS

This cable allows connection of older Bovie electrosurgical units to bipolar forceps.

This cable allows connection of older Bovie electrosurgical units to bipolar forceps. It is 3M long and suited to devices such as the A900 series.

**ORDER CODES:**  
**BOVA827BP**

[Click For More Information](#)

## OLSEN 184MM GERALD SMOOTH 1.0MM TIPS INSULATED BIPOLAR FORCEP

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Olsen 184mm Gerald Smooth 1.0mm are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### Features:

High quality materials  
Insulated  
Suitable for use with all electrosurgical devices

**ORDER CODES:**  
**BOV10-11801**

[Click For More Information](#)



## **OLSEN 140MM SEMKEN SMOOTH 0.75MM TIPS INSULATED BIPOLAR FORCEP**

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Olsen 140mm semken smooth 0.75mm tips are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### **Features:**

High quality materials  
Insulated  
Suitable for use with all electrosurgical devices

**ORDER CODES:**  
**BOV10-10911**

[Click For More Information](#)

## **SYMMETRY 101MM JEWELERS 0.4MM TIPS INSULATED BIPOLAR FORCEP**

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Symmetry 101mm jewelers 0.4mm tips are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### **Features:**

High quality materials  
Insulated  
Suitable for use with all electrosurgical devices

**ORDER CODES:**  
**BOV08-00011**

[Click For More Information](#)





## **OLSEN 127MM ADSON SERRATED 1.0MM TIPS INSULATED BIPOLAR FORCEP**

Reusable Bipolar Forceps are available in micro .5mm, 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8".

Olsen 127mm adson serrated 1.0mm tips are insulated bipolar forceps manufactured to industry gold standard utilising the highest quality materials.

### **Features:**

High quality materials

Insulated

Suitable for use with all electrosurgical devices

### **ORDER CODES:**

**BOV10-10811**

[Click For More Information](#)

# SMOKE EVACUATORS & ACCESSORIES

## BOVIE SMOKE SHARK II SMOKE EVACUATOR



The Bovie Smoke Shark II is designed for ease of use and effective smoke and particle removal at the operative site. The Smoke Shark II is designed for use with all electrosurgical generators and effectively controls smoke plume from lasers and other forms of advanced energy.

The Bovie Smoke Shark II is designed for ease of use and effective smoke and particle removal at the operative site. The Bovie Smoke Shark II is designed for use with all electrosurgical generators and effectively controls smoke plume from lasers and other forms of advanced energy.

### Features:

- Easy to operate
- Quiet and effective smoke evacuation
- Long life 35 hour filter
- Low cost per procedure
- Excellent particle capture

### Specification:

#### Recommended Accessories

- 786T - Vacuum hose 7/8" x 6"
- 786TS - Vacuum hose 7/8" x 6" hose to hose connect
- SS95 - Smoke evacuation pencil adaptor 3/8" x 20" hose
- SERS2 - Ez link smart remote activator switch

#### Dimensions

- Height: 15.2cm
- Width: 27.9cm
- Depth: 39.4cm
- Weight: 4.3KG

#### Line Voltage

220VAC

#### Line Frequency

50-60Hz

#### Warranty

2 year warranty

#### Kit Inclusions

- SE02-220 smoke shark 2 x 1
- FSSEP pneumatic footswitch x 1
- SF35 long life 35 hour filter x 1
- AC power cord x 1
- User guide x 1

**ORDER CODES:  
BOVSE02-220**


**MAINTENANCE:** Please make sure your device is checked and serviced yearly by an authorised technician.

[Click For More Information](#)

## BOVIE SMOKE SHARK III SMOKE EVACUATOR

The Bovie Smoke Shark III is designed for ease of use and effective smoke and particle removal at the operative site. The Smoke Shark III is designed for use with all electrosurgical generators and effectively controls smoke plume.

Bovie smoke shark III reduces biological, physical, and chemical hazards present in electrosurgical smoke plume. Uniquely designed for the physician's office with compact size for wall mounting. Specifically intended for shorter procedures (15 minute or less) with a 4-Stage ULPA filtration 99.999% efficient in trapping 0.1-0.2 micron particulates and micro-organisms. A quiet design offers excellent performance with less noise. Includes a shark skin adapter (SS95) for out of the box usability, a smart shark remove activation device (SERS2) for maximizing filter life & a 6 hour filter (SF06).

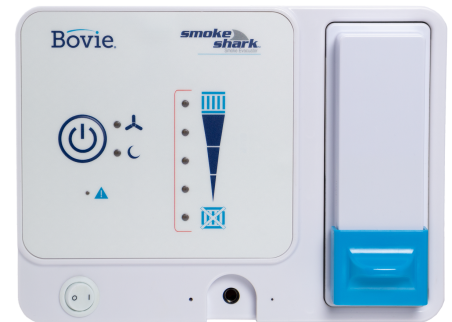

### Features:

- Compact for the physician office
- Quiet and effective smoke evacuation
- Multi procedure 6 hour filters
- Excellent particle capture

### Specification:

#### Dimensions

- Width: 23.9cm
- Depth: 15.5cm
- Height: 18cm
- Weight: 2.8kg

#### Line Voltage

220VAC

#### Line Frequency

50-60Hz

#### Warranty

2 year warranty

#### Kit Inclusions

- SE03-220 smoke shark 3 x 1
- SERS2 ezi link remote activator x 1
- SS95 smoke evacuation pencil adaptor x 1
- SF06 6 hour filter x 1
- Wall mounting kit x 1
- AC power cord x 1

User guide x 1

**ORDER CODES:**  
**BOVSE03-220**

**MAINTENANCE:** Please make sure your device is checked and serviced yearly by an authorised technician.

[Click For More Information](#)



## **BOVIE SMOKE SHARK II 35 HOUR FILTER NON-STERILE**

Bovie smoke shark II 35 hour filter suitable for use with the Smoke Shark II. Designed for ease of use and effective smoke and particle removal at the operative site.

Bovie smoke shark II 35 hour filter suitable for use with the Smoke Shark II. Designed for ease of use and effective smoke and particle removal at the operative site.

### **Features:**

Easy replacement  
35 hour filter life

**ORDER CODES:**  
**BOVSF35**

[Click For More Information](#)

## BOVIE SMOKE SHARK III FILTER NON-STERILE

Bovie smoke shark III 6 hour filter suitable for use with the Smoke Shark III. Designed for ease of use and effective smoke and particle removal at the operative site.

Bovie smoke shark III 6 hour filter suitable for use with the Smoke Shark III. Designed for ease of use and effective smoke and particle removal at the operative site.

### Features:

Easy replacement  
6 hour filter life

**ORDER CODES:**  
**BOVSF06**

[Click For More Information](#)



## BOVIE SMOKE EVACUATOR PNEUMATIC FOOT SWITCH

Smoke Evacuator Pneumatic Footswitch to suit use with smoke shark devices.

Smoke Evacuator Pneumatic Footswitch to suit use with smoke shark devices.

**ORDER CODES:**  
**BOVFSSEP**

[Click For More Information](#)

**BOVIE VACUUM HOSE 7/8" X 6" BX24**

Vacuum hose to hose connect 7/8" x 6" to suit use with smoke shark devices.

Bovie vacuum hose 7/8" x 6" to suit use with smoke shark devices.

**ORDER CODES:**  
**BOV786T**



[Click For More Information](#)

**BOVIE REDUCER FITTER 7/8" TO 1/4" X 24" BX10**

Bovie Reducer Fitter 7/8" to 1/4" x 24" to suit use with smoke shark devices.

Bovie Reducer Fitter 7/8" to 1/4" x 24" to suit use with smoke shark devices.

**ORDER CODES:**  
**BOVSERF**

[Click For More Information](#)

## **BOVIE SMOKE EVACUATION PENCIL ADAPTOR 3/8" X 20" HOSE BX20**

Smoke evacuation pencil adaptor 3/8" X 20" hose to suit use with smoke shark devices.

Smoke evacuation pencil adaptor 3/8" X 20" hose to suit use with smoke shark devices.

**ORDER CODES:  
BOVSS95**



[Click For More Information](#)

## **BOVIE EZ LINK**

Bovie ez link remote activator switch to suit use with smoke shark devices.

Bovie ez link remote activator switch to suit use with smoke shark devices. Simply clip onto the handpiece cable to allow automatic activation of the Smoke Shark smoke evacuator.

**ORDER CODES:  
BOVSERS2**

[Click For More Information](#)



## **BOVIE SMOKE EVACUATOR ARTICULATED ARMSTAND**

Smoke Evacuator Articulated Armstand to suit use with smoke shark devices.

Smoke evacuator articulated armstand to suit use with smoke shark devices.

**ORDER CODES:**  
**BOVFSEAS**



[Click For More Information](#)

## **BOVIE HANDHELD VACUUM HOSE 7/8" X 6" WITH WAND & SPONGE GUARD BX24**

Handheld vacuum hose 7/8" x 6" with wand & sponge guard to suit use with smoke shark devices.

Handheld vacuum hose 7/8" x 6" with wand & sponge guard to suit use with smoke shark devices.

**ORDER CODES:**  
**BOVSETWS**

[Click For More Information](#)



**BOVIE VACUUM HOSE 7/8" X 6" HOSE TO HOSE  
CONNECT BX24**

Vacuum hose to hose connect 7/8" x 6" to suit use with smoke shark devices.

Vacuum hose to hose connect 7/8" x 6" to suit use with smoke shark devices.

**ORDER CODES:**  
**BOV786TS**

[Click For More Information](#)



# MISCELLANEOUS & SPARE PARTS



## BOVIE COBALT BLUE DISPOSABLE PENLIGHT PK3

Cobalt blue disposable penlight are used in conjunction with fluorescein dye to find foreign matter in the eye.

Cobalt blue disposable penlight are used in conjunction with fluorescein dye to find foreign matter in the eye.

**ORDER CODES:**  
**BOV6303**

[Click For More Information](#)

## BOVIE U.V WOODS LIGHT

Bovie U.V woods light are excellent for detecting certain skin disorders.

Bovie U.V woods light are excellent for detecting certain skin disorders.

**ORDER CODES:**  
**BOVUV59**

[Click For More Information](#)



## BOVIE SURCH LITE 381MM STERILE PK3

Bovie surch lite 381mm are for use in procedures where additional light is needed.

Bovie surch lite 381mm are for use in procedures where additional light is needed. Comes sterile and in a pack of 3 units.

**ORDER CODES:**  
**BOVST15**

[Click For More Information](#)

# CAUTERY

## POWER HANDLES & REPLACEMENT TIPS

### BOVIE POWER HANDLE LOW-TEMP



Bovie manufactures more cauteries than anyone in the world. Therefore, Bovie has the ability to put the right cautery in the physician's hand regardless of the procedure. High temperature cauteries are popular for many types of procedures including, pinpoint hemostasis, controlling diffuse bleeders, and many others.

Bovie power handle low-temp replaceable battery and replaceable tip cauteries are the next advancement in battery operated cauteries where reusability is a must. These products are designed for the office, ER, and clinical procedures. The Bovie power handle low-temp can be purchased individually in high or low temperature models with a variety of tips designed to fit your procedural needs.

#### Features:

- Portable & easy to use
- Safe design with recessed power button
- Reusable handle with switch
- Replaceable batteries
- Replaceable tips

#### Specification:

##### Dimensions

Length: 205mm  
Diameter: 19mm

##### Battery

1 x AA

##### Kit Inclusions

HIT0 power handle low-temp x 1  
H100 high-temp fine tip x 1  
AA battery x 1

##### ORDER CODES:

**BOVHIT1**

**MAINTENANCE:** Please make sure your device is checked and serviced yearly by an authorised technician.

[Click For More Information](#)

## BOVIE POWER HANDLE HIGH-TEMP

Bovie manufactures more cauteries than anyone in the world. Therefore, Bovie has the ability to put the right cautery in the physician's hand regardless of the procedure. High temperature cauteries are popular for many types of procedures including, pinpoint hemostasis, controlling diffuse bleeders, and many others.



Bovie power handle high-temp replaceable battery and replaceable tip cauteries are the next advancement in battery operated cauteries where reusability is a must. These products are designed for the office, ER, and clinical procedures. The Bovie power handle high-temp can be purchased individually in high or low temperature models with a variety of tips designed to fit your procedural needs.

### Features:

- Portable & easy to use
- Safe design with recessed power button
- Reusable handle with switch
- Replaceable batteries
- Replaceable tips

### Specification:

#### Dimensions

Length: 205mm  
Diameter: 19mm

#### Battery

2 x AA

#### Kit Inclusions

HIT1 power handle high-temp x 1  
H101 high-temp fine tip x 1  
AA battery x 2

#### ORDER CODES:

**BOVHIT1**

**MAINTENANCE:** Please make sure your device is checked and serviced yearly by an authorised technician.

[Click For More Information](#)

## BOVIE REPLACEMENT TIP LOW-TEMP FINE BX10



Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

### Features:

- Easy to use
- Sterile drape included
- Suitable for power handles

### ORDER CODES: BOVH100

[Click For More Information](#)

## BOVIE REPLACEMENT TIP HIGH-TEMP FINE BX10

Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.



Replacement tip high-temp fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

### Features:

- Easy to use
- Sterile drape included
- Suitable for power handles

### ORDER CODES: BOVH101

[Click For More Information](#)

## BOVIE REPLACEMENT TIP HIGH-TEMP LOOP BX10



Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

Replacement tip high-temp loop are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

### Features:

- Easy to use
- Sterile drape included
- Suitable for power handles

### ORDER CODES: BOVH103

[Click For More Information](#)

## BOVIE REPLACEMENT TIP HIGH-TEMP FINE ELONGATED BX10

Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.


Replacement tip high-temp elongated fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

### Features:


- Easy to use
- Sterile drape included
- Suitable for power handles

### ORDER CODES: BOVH104

[Click For More Information](#)

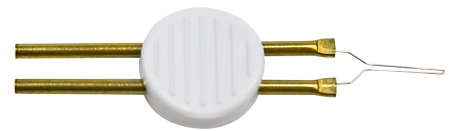

## **BOVIE REPLACEMENT TIP HIGH-TEMP VASECTOMY BX10**

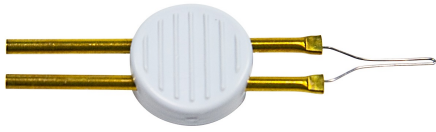

Replacement Tip Low-Temp Fine are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

Replacement tip high-temp vasectomy are sterile cautery tips including a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the power handle cauteries.

### **Features:**

- Easy to use
- Sterile drape included
- Suitable for power handles

### **ORDER CODES: BOVH121**

[Click For More Information](#)

# HIGH TEMPERATURE CAUTERIES

## BOVIE CAUTERY HIGH-TEMP FINE TIP BX10

Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery high-temp fine tip has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.

### Features:

- Portable & easy to use
- Safe design with recessed power button
- Easy snap design

### ORDER CODES: BOVAA01

[Click For More Information](#)



## BOVIE CAUTERY HIGH-TEMP LOOP TIP BX10

Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery high-temp loop tip has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.



**Features:**

Portable & easy to use  
Safe design with recessed power button  
Easy snap design

**ORDER CODES:  
BOVAA03**

[Click For More Information](#)

**BOVIE CAUTERY HIGH-TEMP LOOP TIP EXTENDED  
SHAFT BX10**

Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery high-temp loop tip extended shaft has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.

**Features:**

Portable & easy to use  
Safe design with recessed power button  
Easy snap design

**ORDER CODES:  
BOVAA05**

[Click For More Information](#)

## BOVIE CAUTERY HIGH-TEMP FINE TIP EXTENDED SHAFT BX10

Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.



Bovie cautery high-temp fine tip extended shaft has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.

### Features:

- Portable & easy to use
- Safe design with recessed power button
- Easy snap design

### ORDER CODES: BOVAA17

[Click For More Information](#)

## BOVIE CAUTERY HIGH-TEMP ELONGATED TIP BX10



Bovie has been the number one manufacturer of battery-operated cauterics for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery high-temp elongated tip has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.

### Features:

- Portable & easy to use
- Safe design with recessed power button
- Easy snap design

**ORDER CODES:**  
**BOVAA21**

[Click For More Information](#)

# LOW TEMPERATURE CAUTERIES

## BOVIE CAUTERY LOW-TEMP FINE TIP BX10



Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery low-temp fine tip has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.

### Features:

- Portable & easy to use
- Safe design with recessed power button
- Easy snap design

### ORDER CODES: BOVAA00

[Click For More Information](#)

## BOVIE CAUTERY LOW-TEMP MICRO FINE TIP BX10

Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal.

Bovie cautery low-temp micro fine tip has a safety-enhanced, snap-design & is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to “snap off” the end of the cautery for safe and effective removal of the battery power source.



**Features:**

Portable & easy to use  
Safe design with recessed power button  
Easy snap design

**ORDER CODES:  
BOVAA90**

[Click For More Information](#)

# SLEEVES

## BOVIE CAUTERY STERILE SHEATH BX10

Bovie cautery sterile sheath suitable for use with power handles.



Bovie cautery sterile sheath suitable for use with power handles.

**ORDER CODES:**  
**BOVHISL**

[Click For More Information](#)

# OPHTHALMIC

## BOVIE BURR 1MM BX10

Bovie burr 1mm suitable for use with AB03, AB04 or equivalent.



Bovie burr 1mm suitable for use with AB03, AB04 or equivalent.

**ORDER CODES:**  
**BOVAB01**

[Click For More Information](#)

## BOVIE BURR 1/2MM BX10

Bovie burr 1mm suitable for use with AB03, AB04 or equivalent.

Bovie burr 1/2mm suitable for use with AB03, AB04 or equivalent.

**ORDER CODES:**  
**BOVAB05**

[Click For More Information](#)



# INDEX

## A

Autoclavable Three Button Pencil New 11

Autoclavable Three Button Pencil Old 12

## B

Bovie 114mm Straight Micro Tip Bipolar Forcep 38

Bovie 177mm Straight Fine Tip Bipolar Forcep 38

Bovie 184mm Fine Tip Bipolar Forcep 37

Bovie Autoclavable Reusable Pencil 12

Bovie Bantam PRO High Frequency Desiccator 5

Bovie Burr 1/2mm BX10 65

Bovie Burr 1mm BX10 65

Bovie Cautery High-Temp Elongated Tip BX10 60

Bovie Cautery High-Temp Fine Tip BX10 58

Bovie Cautery High-Temp Fine Tip Extended Shaft BX10 60

Bovie Cautery High-Temp Loop Tip BX10 58

Bovie Cautery High-Temp Loop Tip Extended Shaft BX10 59

Bovie Cautery Low-Temp Fine Tip BX10 62

Bovie Cautery Low-Temp Micro Fine Tip BX10 62

Bovie Cautery Sterile Sheath BX10 64

Bovie Cobalt Blue Disposable Penlight PK3 51

Bovie Derm 102 High Frequency Desiccator 2

Bovie Derm 942 High Frequency Desiccator 3

Bovie Disposable Derm-Elite Round Electrode Tips Blunt Angled Non-Sterile BX100 22

Bovie Disposable Derm-Elite Round Electrode Tips Blunt Angled Sterile BX50 22

Bovie Disposable Dermal Electrode Tips Blunt Angled Non-Sterile BX100 21

Bovie Disposable Dermal Electrode Tips Blunt Angled Sterile BX50 21

Bovie Disposable Dermal Electrode Tips Sharp Non-Sterile BX100 20

Bovie Disposable Dermal Electrode Tips Sharp Sterile BX50 20

Bovie Disposable Extended Blade Electrodes BX25 24

Bovie Disposable H-Type Electrode Tips Blunt Angled Non-Sterile BX100 29

Bovie Disposable H-Type Electrode Tips Blunt Angled Sterile BX50 30

Bovie Disposable H-Type Electrode Tips Sharp Angled Non-Sterile BX100 28

Bovie Disposable H-Type Electrode Tips Sharp Sterile BX50 29

Bovie Disposable Handpiece Sheath Non-Sterile BX100 15

Bovie Disposable Handpiece Sheath Sterile BX25 15

Bovie Disposable IMA / ENT Extended Needle Electrodes BX25 28

Bovie Disposable IMA / ENT Standard Blade Electrodes BX25 27

Bovie Disposable IMA / ENT Standard Needle Electrodes BX25 27

Bovie Disposable Micro Dissection Needle Electrode BX5 25

Bovie Disposable Push Button Pencil With Holster & Blade BX40 11

Bovie Disposable Restick II Coated Blade Electrodes BX12 23

Bovie Disposable Restick II Coated Modified Needle Electrodes BX12 26

Bovie Disposable Restick II Coated Needle Electrodes BX12 25

Bovie Disposable Split Adult Return Electrode BX50 18

Bovie Disposable Split Adult Return Electrode With 3M Cable BX50 17

Bovie Disposable Standard Blade Electrode BX50 23

Bovie Disposable Standard Needle Electrodes BX25 24

Bovie Disposable Super Cut Modified Fine Needle Electrodes BX5 26

Bovie Ez Link 48

Bovie Foot Switch 9

Bovie Handheld Vacuum Hose 7/8" x 6" With Wand & Sponge Guard BX24 49  
 Bovie Mobile Cart For ESU & Smoke Evacuator 9  
 Bovie Mobile Stand 8  
 Bovie Monopolar / Bipolar Foot Switch 10  
 Bovie Power Cable 10  
 Bovie Power Handle High-Temp 54  
 Bovie Power Handle Low-Temp 53  
 Bovie Reducer Fitter 7/8" to 1/4" x 24" BX10 47  
 Bovie Replacement Metal Grounding Plate 18  
 Bovie Replacement Tip High-Temp Fine BX10 55  
 Bovie Replacement Tip High-Temp Fine Elongated BX10 56  
 Bovie Replacement Tip High-Temp Loop BX10 56  
 Bovie Replacement Tip High-Temp Vasectomy BX10 57  
 Bovie Replacement Tip Low-Temp Fine BX10 55  
 Bovie Reusable 1/4" Tungsten Loop Electrode 34  
 Bovie Reusable 3/32" Ball Electrode 36  
 Bovie Reusable 3/8" Tungsten Loop Electrode 35  
 Bovie Reusable 5mm Extended Ball Electrode 32  
 Bovie Reusable Angled Ball Electrode 33  
 Bovie Reusable Angled Needle Electrode 35  
 Bovie Reusable Bipolar Forceps Cord For New Units 37  
 Bovie Reusable Bipolar Forceps Cord For Old Units 40  
 Bovie Reusable Epilation Needle Electrode 31  
 Bovie Reusable Grounding Cable 19  
 Bovie Reusable Metal Grounding Plate With Cable 17  
 Bovie Reusable Standard Blade Electrode 33  
 Bovie Reusable Standard Needle Electrode 32  
 Bovie Reusable Standard Segal Electrode 34  
 Bovie Reusable Straight Needle Electrode 31  
 Bovie Scratch Pads BX40 19  
 Bovie Smoke Evacuation Pencil Adaptor 3/8" X 20" Hose BX20 48  
 Bovie Smoke Evacuator Articulated Armstand 49  
 Bovie Smoke Evacuator Pneumatic Foot Switch 46  
 Bovie Smoke Shark II 35 Hour Filter Non-Sterile 45  
 Bovie Smoke Shark II Smoke Evacuator 43  
 Bovie Smoke Shark III Filter Non-Sterile 46  
 Bovie Smoke Shark III Smoke Evacuator 44  
 Bovie Specialist PRO High Frequency Electrosurgical Generator 6  
 Bovie Suction Coagulator Disposable Footswitch 10Fr 3M Cable PK10 13  
 Bovie Suction Coagulator Disposable Handswitch 12Fr 3M Cable PK10 13  
 Bovie Suction Coagulator Disposable Handswitch 8Fr 3M Cable PK10 14  
 Bovie Surch Lite 381mm Sterile PK3 52  
 Bovie U.V Woods Light 51  
 Bovie Vacuum Hose 7/8" x 6" BX24 47  
 Bovie Vacuum Hose 7/8" x 6" Hose To Hose Connect BX24 50  
 Bovie Wall Mount Kit 8

## O

Olsen 127mm Adson Serrated 1.0mm Tips Insulated Bipolar Forcep 42  
 Olsen 127mm Adson Smooth 1.0mm Tips Insulated Bipolar Forcep 39  
 Olsen 140mm Semken Smooth 0.75mm Tips Insulated Bipolar Forcep 41  
 Olsen 178mm Cushing Smooth Straight 1.5mm Tips Insulated Bipolar Forcep 39  
 Olsen 184mm Gerald Smooth 1.0mm Tips Insulated Bipolar Forcep 40

## S

Symmetry 101mm Jewelers 0.4mm Tips Insulated Bipolar Forcep 41



**Pacific Medical Australia**

**E** [info@pacificmedical.com.au](mailto:info@pacificmedical.com.au)

**P** 1300 076 332

[www.pacificmedical.com.au](http://www.pacificmedical.com.au)